

Crono 30

AMBULATORY INFUSION PUMP



actual size 1:1

Distributed by:



Intra Pump Infusion Systems
920 Minters Chapel Rd, Suite 200
Grapevine, TX 76051
Tel: 866.211.7867 (PUMP)
www.intrapump.com

Manufactured by:
Canè s.r.l. Medical Technology
Via Pavia 105/i
Rivoli (TO) Italy
www.microjet.it

DEVICE DESCRIPTION

Crono 30 is an ambulatory syringe infusion pump intended for the controlled administration of liquids into the patients. This means that the pump can be worn during the infusion.


Infusion pumps using normal commercial syringes are bulky and cumbersome and thus difficult to use in the true sense of "ambulatory". This can have a negative impact on patient compliance and treatment goals.

Aspecial syringe, an integral part of **Crono 30**, provides an efficient reduction of overall pump size. Whether day or night infusion is required, **Crono 30** offers more conveniences and advantages than other bulkier, heavier pump.

TECHNICAL FEATURES

Pump size	80x47x30 mm.
Weight	125 g (battery included).
Battery	3 V type CR 123 A lithium battery (battery life 4/ 6 months).
Syringe	Dedicated with 30 ml capacity and universal male "luer-lock" connector.
Volume administration	Drug volume can be programmed from 1 to 30 ml.
Delivery time	Programmable from 1 to 99 hours, in 15 minutes steps.
Accuracy	+/-2%.
Available occlusion pressure	4 bar +/-1 .
Data storage	Programmable data are automatically stored in the pump's memory; the data are not lost when the battery is removed.
Safety circuits	To control the regular working of the device, alarming in the event of anomaly and displaying error message.
Tamper resistant feature	Allows the pump to be "locked" by disabling the buttons.
Standard equipment supplied	Fabric holder, elastic belt, collar strap, battery tool, two batteries, user guide.
Manufacturer's guarantee	Canè s.r.l. guarantees that this product is free from defects in materials and workmanship for a period of 2 (two) years beginning from the date of purchase

SYRINGE PREPARATION

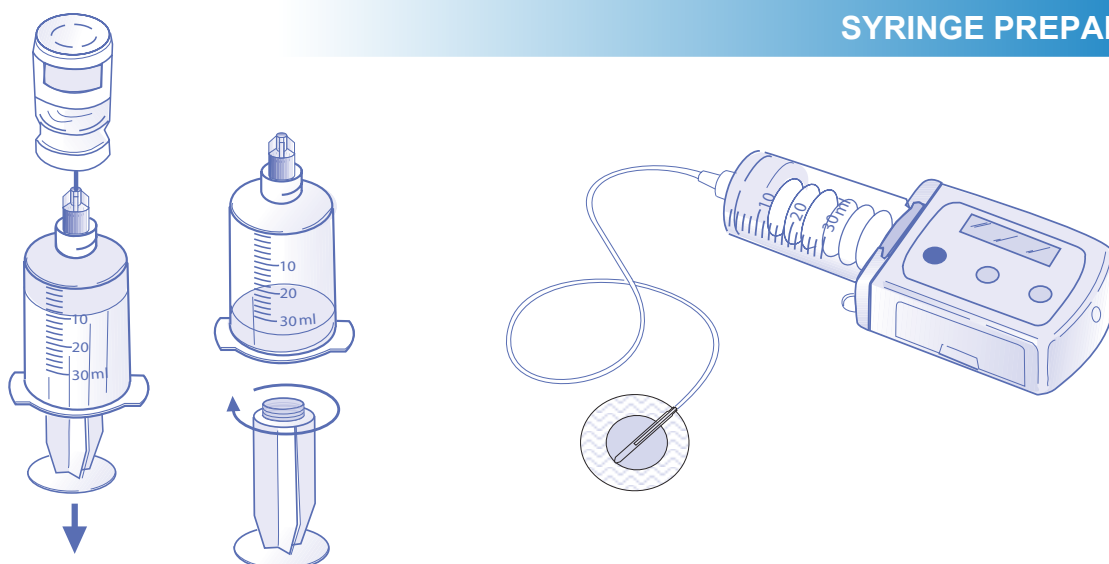

CE
0476